



Primo R

2 in 1 HEARING & MUSIC
RIC HEARING AID




Binaural hearing enhancer



Rechargeable



Music & Phone For All Platform

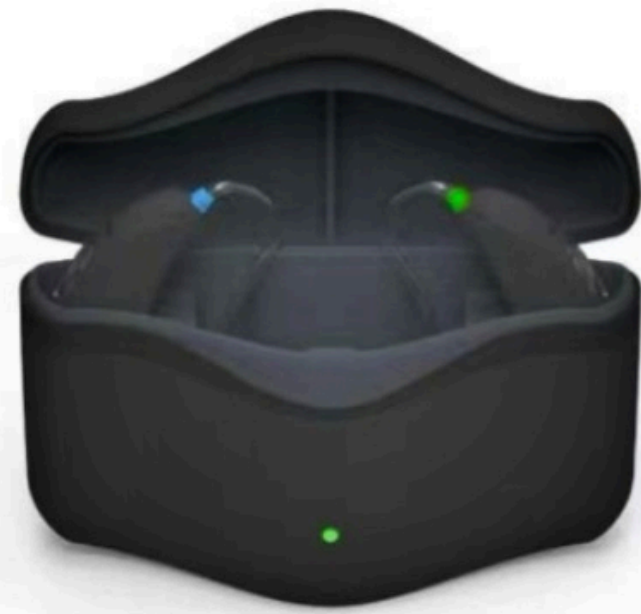

SELLING POINT

2-in-1 Hearing & Music RIC device

Immerse yourself with our comfortable slim BTE design to music streaming and hearing enhancement.

App self-tuning and remote-fitting

Connect wirelessly to iOS® or Android™ devices for SOUNDWEAR® APP remote control: Hearing check / Self tuning Sound Equalizer / Listening Mode Select

FEATURE

- EZCharge® : Magnetic contact charging
- Up to 32 channels WDRC
- Long use time
- Low battery warning
- Drop and Charge, Pick and On
- Intelligent noise reduction function
- User friendly operation button
- 4 listening programs
- Resume of memory
- Dual microphones for better speech understanding

PACKAGE INFO

- Type-C charging cable
- Silicone ear dome
- Charging case
- Brush

TECHNICAL DATA

Maximum output OSPL 90(dBSPL)	120
HFA output OSPL 90(dBSPL)	114
Maximum Full-on gain(dB)	52
HFA Full-on gain(dB)	43
THD(%)	< 2
EQ Input Noise(dB)	< 28
Frequency Range (Hz)	300-5000
Battery Current (mA)	4.5

Battery type	Lithium
Working hours (Hr) Hearing	10 hours
Working hours (Hr) Music	7 hours
Charging time (per pair)	3 times
Time of full charge (Devices)	2.5 hours
Time of full charge (Base)	3 hours
Size and weight (Device)	32.9×13×8.9mm
Size and weight (Base)	80.1×54.7×32mm

Data in accordance with ANSI S3.22-2003 IEC 60118-7:2005 2cc coupler.
All technical(typical) data are subject to change without prior notice.

IMHear Corporation

Add: 2711 Curtiss Street, Suite B, Downers Grove, IL, 60515, USA Tel: 1-630-395-9628
Website: www.imhearing.com Email: info@imhear.comcastbiz.net